














Chattanooga Primera

DUAL CHANNEL TENS & NMES UNIT

Operators Manual



Symbols on the unit and case

	Caution! (electrical output).
	Follow operating instructions! Failure to do so could place the patient or operator at risk.
	Neuromuscular Stimulation (STIM) and EMG Triggered Stimulation (ETS) should not be used by Patients fitted with demand style cardiac pacemakers. Please seek advice from your health supervisor.
	Patient's shock protection type: BF (Body floated) Equipment. Floating isolated applied part. It is only intended for connection to patient's skin but has floating input circuits. No connections between patient and earth.
	Indicates the manufacturer's catalogue number so that the medical device can be identified.
	Manufacturer's LOT/Batch number. Present it together with SN number when you report a technical fault or claim a warranty return.
	Manufacturer's serial number of the unit. Present it together with LOT number when you report a technical fault or claim a warranty return.
	Name and address of Manufacturer.
	Date of manufacture.
	Conformity indication with the essential health and safety requirements set out in European Directives. 0086 - Notified Body identification (BSI).
	The Australian government requires that all imported or locally produced electrical and electronic equipment comply with electromagnetic compatibility (EMC) emission requirements. A compliant product must bear the C-Tick logo.
	This product should be kept dry.
IP20 on the unit	This is an indication for protection against ingress of water and particulate matter. The mark IP20 on your unit means: your unit is protected against solid foreign objects of 12.5mm dia and greater. Not protected against water.
IP02 on the case	IP02 on the carrying case means: Protected from the ingress of water droplets from a shower of rain.
	Do not dispose in normal dustbin (see page 20 for the disposal instructions).

Warnings



- * This unit must be used with the guidance of a Physiotherapist or Doctor.
- * Type BF equipment, Continuous Operation.
- * Do not insert lead wires into a mains power supply.
- * Do not immerse unit into water or any other substance.
- * The unit is not protected from the ingress of water droplets from a shower of rain if used outside the carrying case.
- * Do not use this unit in the presence of a flammable anaesthetic gas mixture and air or with Oxygen or Nitrous Oxide.
- * If using rechargeable 9 Volt PP3 Nickel Metal Hydride batteries, be sure to use a CE approved battery charger. Never connect this unit directly to a battery charger or to any other mains powered equipment.
Do not use Ni-Cad rechargeable batteries.
Caution: Do not use lithium batteries unless they comply with IEC60086-4.
- * Patient Electrodes are for single patient use only.
- * Keep out of reach of children.
- * Do not use this stimulator on your facial area unless you are under strict guidance from a qualified Clinician.
- * Application of electrodes near the thorax may increase the risk of cardiac fibrillation.
- * Operation in close proximity (e.g. 1m) to a shortwave or micro wave therapy equipment may produce instability in the stimulator output.
- * Simultaneous connection of a patient to a high frequency surgical equipment may result in burns at the site of the stimulator electrodes and possible damage to the stimulator.
- * No modification of this equipment is allowed!



Contents

Contents	Page
Symbols of the unit	2
Warnings	3
Intended purpose	5
What is Pain?	5
What is TENS?	5
What is STIM?	6
Contraindications & Precautions	7
Description of Unit & Functions	8
Quick Start Instructions & Important Notes	9
Programmes	10
Using the Chattanooga Primera™ Unit in TENS mode	11
Treatment Modes	12
Electrode Placement (TENS)	13
Electrode Types & Tips	14
Electrode Positioning (TENS)	15
Electrode Positioning (NMES)	16
Care, Maintenance, Accessories and Disposal	20
Accessories and re-order codes	22
Conditions that Respond to TENS	23
Conditions that Respond to STIM (NMES)	23
Information regarding electromagnetic compatibility and interference (EMC)	24
Specifications	28
Issues and Solutions	29
Warranty	32
Dermatomes & Myotomes	33
Clinical References	35



Intended Use

Chattanooga Primera™ uses a small battery operated unit to provide a non-invasive, drug free method of controlling acute and principally long term intractable pain. It can also be used as adjunctive treatment in the management of post surgical traumatic pain problems. Mild electrical impulses are transmitted through the skin via surface electrodes to modify the body's pain perception and/or to stimulate muscle and nerve fibres using precise electrical signals via various predetermined or programmable stimulation modulation options.

What is Pain?

When we feel pain it is the body's process of informing us that something is wrong. To feel pain is important, without this feeling abnormal conditions may go undetected, creating damage or injury to critical parts of the body.

Although pain is essential in warning our body of trauma or malfunction, nature may have gone too far in its design. Continued long-term chronic pain has no useful value apart from its importance in diagnosis. Pain begins when a coded signal travels to the brain where it is decoded, and analysed. The pain message travels from the injured area of the body along small diameter nerves leading to the spinal cord. At this point the message is switched to a different kind of nerve that travels up the spinal cord to the brain area. The brain then analyses the pain message, refers it back and the pain is felt.

What is TENS?

Transcutaneous Electrical Nerve Stimulation (TENS) uses a small battery operated unit to provide a non-invasive, drug free method of controlling acute and principally long term intractable pain. It can also be used as an adjunctive treatment in the management of post surgical traumatic pain problems. In TENS mild electrical impulses are transmitted through the skin via surface electrodes to modify the body's pain perception. TENS does not cure problematic physiological conditions; it only helps to control the pain perception. TENS will not work for every user, please seek advice from your Doctor.

There are millions of small nerve fibres throughout the body and it only requires a few impulses to produce chronic pain. In addition to small fibres, which allow the sensation of pain to be felt, the body is also made up of larger diameter nerve fibres. These larger nerve fibres transmit less unpleasant sensations such as touch or warmth, assisting us to form an impression of our environment. Stimulating the larger nerve fibres using TENS may have the effect of inhibiting the transmission of pain along the smaller nerve fibres to the spinal cord [known as the 'Pain Gate Theory'].



What is STIM?

Neuromuscular Stimulation has been used for many years to stimulate muscle and nerve fibres to treat a number of muscle and nerve related conditions. Over the last 30 years numerous clinical trials and papers have been written.

Chattanooga Primera™ is a dual channel device combining several treatment programmes into one unit. Neuromuscular Stimulation is increasingly understood by Therapists and Doctors. There is a better understanding of the mechanisms which exist between nerves and muscles that makes it possible to stimulate the neuromuscular system with precise electrical signals.

Chattanooga Primera™ offers precision giving full control of Pulse Widths, Rates, Ramp up times, Work / Rest cycles as well as alternating or synchronous application if two channels are being applied.

Customer Care

We welcome constructive comments regarding our equipment particularly those that might help us to improve existing features, add new ones or develop new products for the future.

Contraindications & Precautions




Before using this equipment you must first seek the advice of your Physiotherapist or Doctor.

Read this operating manual before using the unit


TENS and STIM should not be used:

- * By patients fitted with a demand style cardiac pacemakers unless so advised by their Doctor.
- * During pregnancy [unless medically advised].
- * By patients with undiagnosed pain conditions.
- * By patients with undiagnosed skin conditions.
- * With patients who have diminished mental capacity or physical competence who cannot handle the device properly.
- * On anaesthetised or desensitised skin.
- * When driving a vehicle or operating potentially dangerous equipment.
- * Do not place electrodes:
 - > Over carotid sinus nerves
 - > Over larynx or trachea.
 - > Inside mouth.
 - > Over the area of the heart unless so advised by your Doctor.
 - > On your facial area unless under strict guidance from a qualified Clinician.
 - > Do not apply stimulation across or through the head, directly on the eyes, covering the mouth, on the front of the neck (especially the carotid sinus) or via electrodes placed on the chest and upper back or crossing over the heart.
- * The patient should use the unit only as prescribed.
- * Do not immerse the unit in water or any other liquid.
- * If you experience skin irritation this may be due to over-stimulation. In this case leave the skin to heal and use TENS only for the periods prescribed. Turning the current up too high can cause skin irritation. In this case allow the skin to heal and use TENS at a lower intensity. Some people experience an allergic reaction to the adhesive coating on the surface of the electrode. If this happens use a different make of electrode or change the electrode. If it continues try reducing the pulse width. If the problem still persists try moving the electrode position each day by just the width of the electrode, making sure the electrode positioning is still over the dermatome.
- * Keep unit out of reach of children.
- * Only use CE approved skin electrodes.
- * If in doubt about the use of the Primera unit, call your Doctor, Therapist, Clinician or you distributor for advice.

Description of Unit & Functions

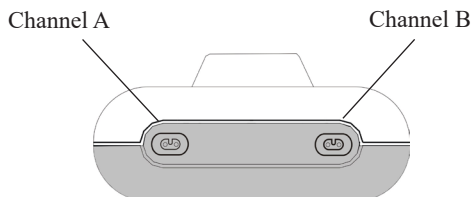

Front

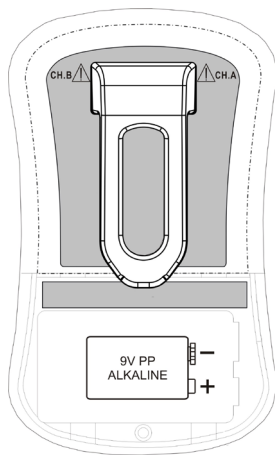

Back

*** PRG button**

Press the PRG [Programme] button to select:
 P01 - P07 or HAN for preset TENS programmes.
 P08 - P13 for preset STIM programmes.
 The summary of all programmes is on page 10.

*** ON / OFF button**

Turns unit on, off and ends the current programme



Quick Start Instructions

1. Insert a 9 volt PP3 Alkaline battery. Alternatively insert a rechargeable Nickel Hydride battery. Do not use Ni-Cad rechargeable batteries
2. Insert lead wire/s to channel A and B if both channels are to be used.
3. Switch on the unit by pressing the ON/OFF button.
4. Press the PRG [Programme] button to select:
P01 - P07 or HAN for preset TENS programmes.
P08 - P13 for preset STIM programmes.
The summary of all programmes is on page 10.
5. To start press channel A + and B + button if you are using both channels, increase the stimulation to the desired level.
6. To stop the programme, press the ON/OFF button which will turn the unit off.

Important Notes

Please carefully read this User's Manual before using your **PRIMERA™**.

Electrode must be in position before the **PRIMERA™** is turned "ON."

Electrodes are for single patient use only - do not share electrodes with another person.

Only use a 9Volt (standard 800 mAh) Alkaline battery in your **PRIMERA™** device.

Do not use Ni-Cad rechargeable batteries or Lithium batteries in the **PRIMERA™** device.

The **PRIMERA™** is equipped with a milliamp (mA) lockout. Forty-five (45) seconds after the **PRIMERA™** is turned "ON," the mA level of stimulation will lock at the user's last mA setting. This prevents an inadvertent increase in the mA setting during use. To increase the mA level forty-six (46) seconds or more after turning the **PRIMERA™** "ON," you must press the negative (-) button first and lower the mA level two (2) mA before increasing it.

The lead wires fit snugly into the controlling unit receptacles for Channels A & B. Please do not twist the lead wires plug while inserting or removing the lead wire plug because twisting may damage the plug (damage may not always be visible). Replace the lead wires regularly to maintain optimum performance of your device.

Programmes

	Programme Number	Programme Description	Rate (Hz)	Pulse width (μ S)	Work time (s)	Rest time (s)	Programme time
TENS	P01	CON	80	200	n/a	n/a	1 hour
	P02	CON	150	175	n/a	n/a	20 min
	P03	CON	2	175	n/a	n/a	20 min
	P04	BST	150	200	n/a	n/a	1 hour
	P05	BST	150	175	n/a	n/a	1 hour
	P06	MOD	100/65	200/100	n/a	n/a	1 hour
	P07	MOD	65/100	200/100	-	-	1 hour
	HAN	MOD	2/70	250/150	-	-	30 min
NMES	P08	NMES	12	200	5	5	15 min
	P09	NMES	35	250	8	8	15 min
	P10	NMES	12	200	5	10	15 min
	P11	NMES	35	200	6	12	15 min
	P12	NMES	12	250	5	15	15 min
	P13	NMES	35	200	6	18	15 min

Key to Program Description for above

CON = CONTINUOUS TENS

BST = BURST TENS

MOD = MODULATION TENS

HAN = HAN TYPE OF MODULATION TENS

NMES = NEURO MUSCULAR ELECTROSTIMULATION



Using the Chattanooga Primera™ unit in TENS mode

RATE [Hz or pulses per second]

The **RATE** to be selected depends primarily on the electrode placement on the patient's body. If one uses contiguous and dermatome (the electrodes alongside or over the area of pain) electrode placement, a higher rate of 80 Hz - 100 Hz is desirable. The patient should experience steady continuous stimulation. It has been found that an optimal setting of 80 or 90 Hz with a pulse width of 200 μ S has good effect for most patients and is a good first choice for pain-gating. Patients using Trigger, motor or acupuncture points tend to respond to low rate stimulation 2 Hz - 10 Hz and pulse width of 200 μ S. The desired effect is for the patient to feel individual pulses.

PULSE WIDTH [Duration]

The wider pulse widths will deliver stronger stimulation for any given intensity [mA] setting. By using a combination of intensity and pulse duration, it is felt that various pulse widths are capable of stimulating different groups of nerve fibres. The wider pulse duration is needed to recruit motor fibres, where as the narrow pulse duration is used more on the sensory fibres. The selection of which pulse duration to use is dependent upon the intended treatment protocol.

Stimulating the larger nerve fibres is thought to reduce the speed and the amount at which information is transmitted along the smaller nerve fibres. Also under certain circumstances the brain is thought to produce its own analgesic pain-killing substances, known as endorphins or endogenous opioids.

Intensity [mA]

Patients respond differently to the level of intensity, this is due to differences in individual patient's skin resistance, enervation and the type and condition of electrode being used.

A good formula for setting the intensity is to increase the current so that the patient feels slight muscle contraction, but not strong enough to move a joint, and then slightly reduce the intensity so that it feels comfortable. When using low rate TENS settings, individual twitches will occur. The higher rate TENS settings will increase muscle tension. It is not advised to increase the intensity to experience strong muscle contraction.



Treatment Modes (TENS)

There are four treatment modes available on the Chattanooga Primera™ unit:

- 1. Conventional TENS or normal.** Train of constant value pulses on both channels delivered with no interruption for entire treatment period. This is the most frequently used mode of the four TENS modes. The most common selection is 80 Hz with a 200 μ S pulse width.
- 2. Burst Mode.** This mode is comparable to the low rate TENS technique except that each low rate pulse is substituted for by a short BURST of 9 pulses [200 μ S] at 150 Hz. It is a combination of conventional and low rate TENS. The burst mode is often referred to as acupuncture - like TENS.
- 3. Modulation TENS.** This mode was designed to help prevent nerve accommodation that some patient's experience. It is achieved by continuously cycling the pulse width and rate.
- 4. Han TENS.** Han TENS is a variation of modulated TENS. It is described as sequence of low and high frequency of stimulation [Dense-and-Disperse (DD) modes of stimulation] where 2 Hz is alternating with 70 Hz, each period lasting for 3 seconds.

How Long Do I Use TENS For?

This depends on the individual patient's condition, accuracy of electrode placement, stimulation and the characteristics selected, but typically the onset of pain relief starts after 20 - 30 minutes. Generally TENS is used for longer periods of normally 1 hour 30 minutes per session. With some patients it can be much longer.



Electrode Placement (TENS)

The placement of electrodes is one of the most important parameters in achieving effective pain relief using TENS. This is best left to your Physiotherapist or Doctor to advise as to which location is most appropriate. It may transpire that various positions need to be experimented with before the user finds the most effective positioning. The positioning may be via the contiguous, dermatome, myotome, motor, trigger or acupuncture points.

Dermatomes & Myotomes

These are areas of the body enervated by a single nerve root via the spinal cord. Each nerve root serves a known area of the skin. The dermatomes are named after the nerve root which serves it. For details of dermatome sites refer to diagrams on pages 33 & 34.

Contiguous Placement

This form of electrode placement is the most common method used. It involves placing the red lead [proximal] alongside the spine where the dermatome [on which your pain lies] enters and exits. The black lead [distal] is normally placed over or near to the pain site. Your Physiotherapist or Doctor may direct the current to cross through the pain area or using the 'bracket' system allow the current to flow on either side of the pain site through the nerve branches that supply the pain location.

Acupuncture Points

The placement of the red and black electrodes on the skin forms the electrical circuit for TENS. It is the skin itself that creates the highest electrical resistance to stimulation. The Physiotherapist or Doctor may consider using acupuncture loci, which offer much lower resistance properties, as a more effective site for placing the electrodes.

Accurately locating an acupuncture point can be difficult, please seek advice from you doctor or physiotherapist.




Electrodes Types and Tips

- * Self-Adhesive reusable long-term electrodes (with proper care) have a typical life span of 4/6 weeks. We recommend cleaning the skin with an alcohol-based wipe before placing the electrodes. The wipe should not contain fat as any grease will degrade the electrode stickiness. After use, place the electrodes back onto the plastic film and in the zip-tag plastic pouch. Store in a cool environment which is not too dry.

A Few Good Tips [Self- Adhesive Electrodes]

- * If you find the electrodes will not stick due to oily skin, cleanse the skin with soap and water, then rinse and dry the area around the electrode site. If this does not work, try cleansing the skin with a swab impregnated with alcohol.
- * Clip away hairy skin using scissors; don't use a razor to remove the hairs!
- * The electrodes conductive material is water- based. At some point the electrodes will become dry. Moisten the adhesive surface with a few drops of water, and apply onto the plastic film overnight. This procedure will increase the electrode life by few more days.


Electrode Positioning (TENS)

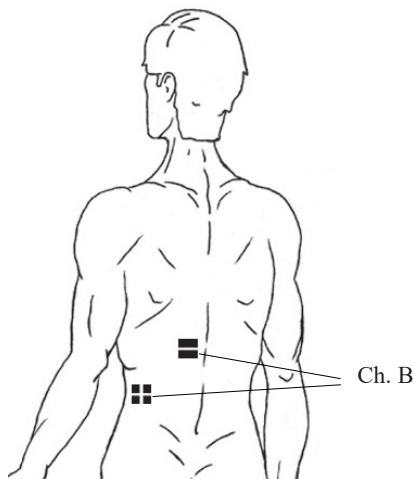

Herpes Zoster

Red = (+) Positive Lead


Black = (-) Negative Lead

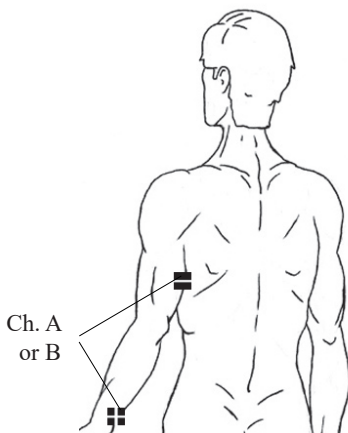

Phantom Limb Pain


Red = (+) Positive Lead

Black = (-) Negative Lead

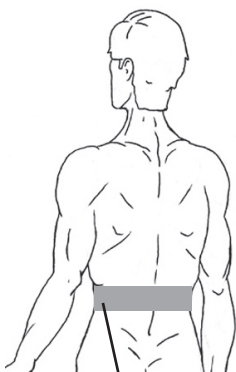


In treating lower-back pain, electrode placement is unique to the user and location should be determined by user response using systematic trial and error. Use of the garment belt (included) is highly recommended to aid in the placement of electrodes on the lower back.

Attach the electrodes to the garment belt and remove the plastic backing. Insert the lead wire on the belt into Channel A or B of the unit. Attach the belt to your body with the electrodes placed on the lower back and secure the belt comfortably around your front. Turn on the unit and increase the current intensity until you feel the desired sensation.

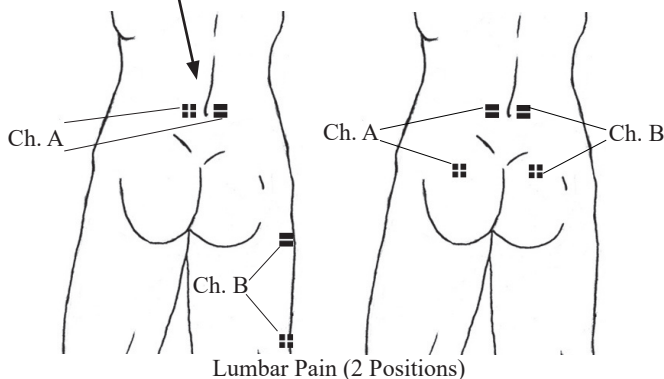


Red = (+) Positive Lead
 Black = (-) Negative Lead

Electrode Positioning (NMES)

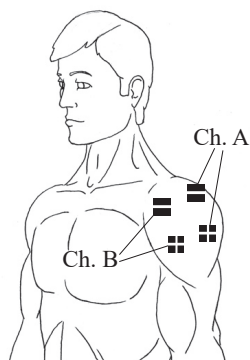


Deltoids

Red = (+) Positive Lead
Black = (-) Negative Lead


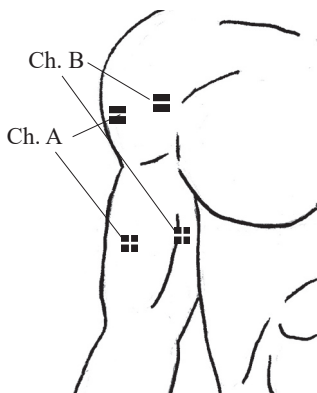


Triceps

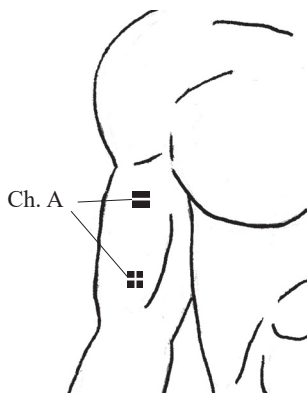


Biceps



Red = (+) Positive Lead
Black = (-) Negative Lead

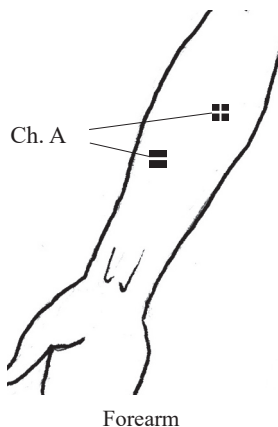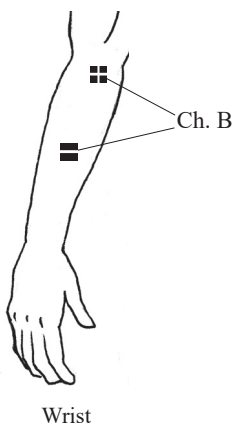



Red = (+) Positive Lead
 Black = (-) Negative Lead

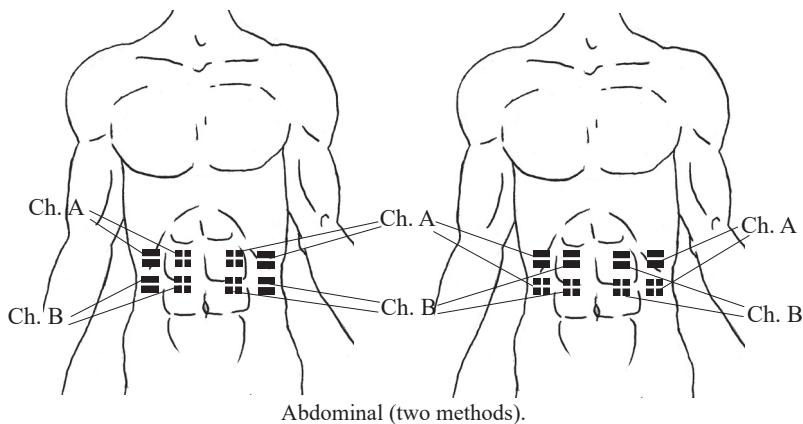

Red = (+) Positive Lead
 Black = (-) Negative Lead

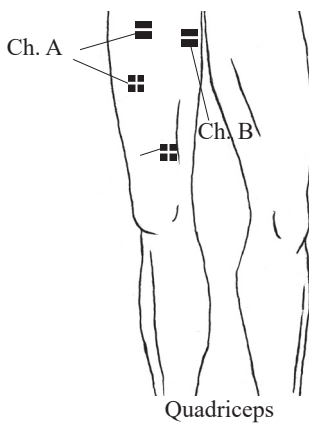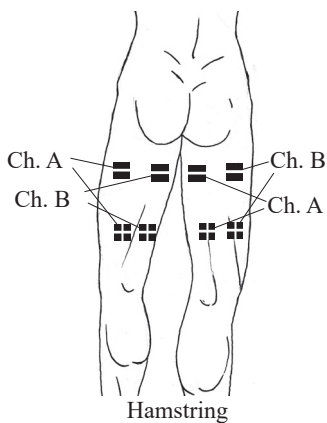

*Red = (+) Positive Lead
Black = (-) Negative Lead*

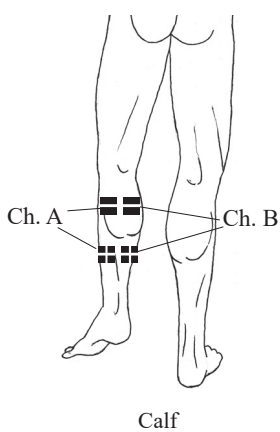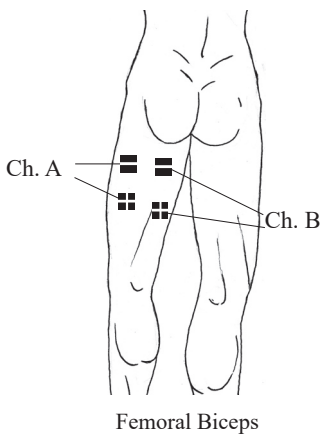

*Red = (+) Positive Lead
Black = (-) Negative Lead*



Care, Maintenance, Accessories and Disposal

WARNING! Only accessories supplied by DJO should be used!

CONTROL UNIT

- * Wipe the surface after each use with a damp cloth or antiseptic wipe or baby wipe.
- * Do not use cleaning sprays or alcohol based cleaning solutions.
- * Expected service life is five years.
- * Control unit disposal: must be disposed of in compliance with national regulatory requirements.

ACCESSORIES

Battery:

- * To change the battery, open the battery door on the rear of the control unit by pressing down on the raised rib pattern just below the belt clip. Lift the battery out of the compartment. This is very easy and can be done by the user.
- * Check periodically for any discharge from the battery.
- * Remove battery completely from unit if not in use for any extended period of time (typically one week).
- * Low battery indicator of 6.9 volts shown on LCD display, when flashing change battery for a new one.
- * Preferably use a PP3 alkaline battery.
- * Expected service life [of a standard 800 mAh, alkaline battery] is 24 hours.
- * Battery disposal: must be disposed of in compliance with national regulatory requirements.

Lead Wires:

- * The lead wires should be handled carefully and never stretched, as this can cause the stimulation to function below normal standards or not at all.
- * Examine lead wires before each treatment for loose connections or damage.
- * Avoid stretching and twisting the lead wires.
- * Store the lead wires carefully after each use.
- * Expected service life is six months, if looked after carefully.
- * Lead wires Disposal: must be disposed of in compliance with national regulatory requirements.



Self-Adhesive Electrodes:

- * Check the short connectors have not become separated from the electrodes.
- * Replace electrodes onto plastic film after use. If they drop onto the floor debris will adhere to conductive gel rendering the electrodes ineffective.
- * Expected service life is a ten day period within two years of the manufacture date. Expiration date is clearly marked on each electrode package.

Electrode life can be considerably reduced by:

- * The type and condition of the skin.
- * Deep seated moisturisers or make-up.

User's Manual:

- * Keep the user's manual clean, dry and away from any open flame or heat source.

For the Best Results:

- * Before each use cleanse the skin.
- * After each use stick the pads on the shiny insert card and store in a cool and dry place, such as the fridge (not freezer).

Caution: Static electricity may damage this product

***NOTE:** Only DJO or appointed distributors / importers are approved to undertake servicing.*



Accessories and Re-order codes

Description	Item code
Chattanooga Primera™ TENS/NMES Complete Unit (International Contains: Main unit, Lead wires (2), Battery, Electrodes in the Plastic Case, Electrode Positioner Belt for the Back, English IFU (Spanish, French, German and Italian languages are for online download in PDF format).	77621

Accessories: You can obtain replacement accessories from your authorized Chattanooga Primera™ dealer:

Description	Item code
Chattanooga Primera™ User's Manual (English)	77622
Personal Carrying Case	77631
Chattanooga Primera™ Unit Sliding Back Cover	77616
Leadwire Set - 2 each	77619
Electrode Positioner Belt for the Back	77620
Battery, Energizer 9-volt	200001-01
Round or Square 5cm (2") Electrodes to be ordered from DJO Global	PN 42192



Conditions that respond to TENS

- * Pain associated with major amputation.
- * Post operative pain.
- * Pain after cardiac surgery.
- * Back pain.

Conditions that respond to STIM (NMES)

- * Upper limb rehabilitation in stroke.
- * Arm function recovery in stroke.
- * Improvement of local blood circulation.

Also used for non medical purposes to:

- * Warm up prior to exercise.
- * Maintain and improve movement.

Information regarding electromagnetic compatibility and interference (EMC)

Our products are designed to produce very low levels of radio frequency (RF) emissions (interference), to be immune from effects of interference produced by other equipment operating in their vicinity and damage due to electrostatic discharge all when operating in a typical domestic and or clinical environment. They are certified to meet the international EMC standard EN60601-1-2. For more information please refer to the tables 201, 202, 204 and 206 overleaf.

Table 201: Guidance and manufacturer's declaration - electromagnetic emission		
This product is intended for use in the electromagnetic environment specified below. The customer or the user of the This product should ensure that it is used in such environment.		
Emission test	Compliance	Electromagnetic environment guidance
RF emission CISPR 11	Group 1	This product uses RF energy only for its internal function. Therefore, its RF emissions are very low and are not likely to cause any interference in nearby electronic equipment
RF emission CISPR 11	Class B	This product is suitable for use in all establishments , including domestic establishments and those directly connected to the public low voltage power supply network that supplies buildings used for domestic purposes
Harmonic emissions IEC 61000-3-2	Not applicable	
Voltage fluctuations/ flicker emissions IEC 61000-3-3	Not applicable	


Table 202: Guidance and manufacturer's declaration
– electromagnetic immunity

This product is intended for use in the electromagnetic environment specified below. The customer or the user of this product should assure that it is used in such an environment, and that precautions regarding that environment are heeded.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Electrostatic discharge (ESD) IEC 61000-4-2	6 kV contact ±8 kV air	±6 kV contact ±8 kV air	Floors should be wood, concrete or ceramic tile. If floors are covered with synthetic material, the relative humidity should be at least 30 %.
Power frequency (50/60 Hz) magnetic field IEC 61000-4-8	3 A/m	3 A/m	Power frequency magnetic fields should be at levels characteristic of a typical location in a typical commercial or hospital environment.

Table 204: Guidance and manufacturer's declaration
– electromagnetic immunity

This product is intended for use in the electromagnetic environment specified below. The customer or the user of this product should assure that it is used in such an environment.

Immunity test	IEC 60601 test level	Compliance level	Electromagnetic environment guidance
Conducted RF IEC 61000-4-6	3 Vrms 150 kHz to 80 MHz	3 Vrms 150 kHz to 80 MHz	<p>Portable and mobile RF communications equipment should be used no closer to any part of this product, including cables, than the recommended separation distance calculated from the equation applicable to the frequency of the transmitter.</p> <p>Recommended separation distance</p> <p>$d = 1.2 \sqrt{P}$ 150 kHz to 80 MHz, $d = 1.2 \sqrt{P}$ 80 MHz to 800 MHz $d = 2.3 \sqrt{P}$ 800 MHz to 2.5GHz</p>
Radiated RF IEC 61000-4-3	3 V/m 80 MHz to 2,5 GHz	3 V/m 80 MHz to 2,5 GHz	<p>where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer and d is the recommended separation distance in meters (m).</p> <p>Field strengths from fixed RF transmitters, as determined by an electromagnetic site survey^a, should be less than the compliance level in each frequency range^b.</p> <p>Interference may occur in the vicinity of equipment marked with the following symbol: </p>

NOTE 1 At 80 MHz and 800 MHz, the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

^a Field strengths from fixed transmitters, such as base stations for radio (cellular/cordless) telephones and land mobile radios, amateur radio, AM and FM radio broadcast and TV broadcast cannot be predicted theoretically with accuracy. To assess the electromagnetic environment due to fixed RF transmitters, an electromagnetic site survey should be considered. If the measured field strength in the location in which this product is used exceeds the applicable RF compliance level above, this product should be observed to verify normal operation. If abnormal performance is observed, additional measures may be necessary, such as reorienting or relocating this product.

^b Over the frequency range 150 kHz to 80 MHz, field strengths should be less than 3 V/m.

Table 206: Recommended separation distances between portable and mobile RF communications equipment and this product

This product is intended for use in an electromagnetic environment in which radiated RF disturbances are controlled. The customer or the user of this product can help prevent electromagnetic interference by maintaining a minimum distance between portable and mobile RF communications equipment (transmitters) and this product as recommended below, according to the maximum output power of the communications equipment

Rated maximum output power of transmitter W	Separation distance according to frequency of transmitter		
	150 kHz to 80 MHz $d = 1.2 \sqrt{P}$	80 MHz to 800 MHz $d = \sqrt{1.2 P}$	800 MHz to 2,5 GHz $d = \sqrt{2.3 P}$
0,01	0.12	0.12	0.23
0,1	0.38	0.38	0.73
1	1.2	1.2	2.3
10	3.8	3.8	7.3
100	12	12	23

For transmitters rated at a maximum output power not listed above, the recommended separation distance d in meters (m) can be estimated using the equation applicable to the frequency of the transmitter, where P is the maximum output power rating of the transmitter in watts (W) according to the transmitter manufacturer.

NOTE 1 At 80 MHz and 800 MHz, the separation distance for the higher frequency range applies.

NOTE 2 These guidelines may not apply in all situations. Electromagnetic propagation is affected by absorption and reflection from structures, objects and people.

Specifications

TENS and STIM(NMES)

1. Dual channel: individually isolated circuits.
2. Amplitude: 0-80 mA into 500 Ohm load ; indication only. Actual mA will tend to be less than indicated due to electrode impedance: at 1000 Ohms load (Electrodes in poor condition) the maximum will be limited to 70 mA, at 1500 Ohms load the maximum will be limited to 65 mA.
3. Type: Constant Current, maximum output voltage 180 Volts +10 / -30 Volts.
4. Waveform: Asymmetrical, rectangular bi-phasic with zero DC current.
5. Selectable pulse width: 100 μ S - 250 μ S [10% accuracy].
6. Pulse Rate selection: in the continuous mode 2 – 150 Hz [5% accuracy].
7. Mode: Continuous, Burst , Modulated or HAN-Modulated.
8. Burst mode: Bursts of 9 pulses [175 μ S or 200 μ S] at 150 Hz, over 2 seconds.
9. Modulation mode: 6-second cycle of concurrent width modulation and pulse repetition rate modulation.
TENS 6: Width starting at 200 μ S and decreasing exponentially to 100 μ S in three seconds and then returning back to 200 μ S in the next three seconds. Rate starting 100 Hz, decreasing exponentially to 65 Hz and then returning to 100 Hz.
TENS 7: Width starting at 200 μ S and decreasing exponentially to 100 μ S in three seconds and then returning back to 200 μ S in the next three seconds. Rate starting 65 Hz, increasing exponentially to 100 Hz over 3 seconds and then returning to 62 Hz in the next 3 seconds.
10. Han mode - variation of Modulated: a sequence of alternating Low frequency stimulation (2Hz, 260 μ S) and High frequency stimulation (70Hz, 150 μ S) each lasting 3 seconds.
11. Battery: PP3 Alkaline, 9V.
Expected average battery life [of standard 800 mAh, alkaline]: 24 h.
12. Low Battery Indicator: If the battery goes below 6.9 volts +/- 0.2 volts the battery symbol will flash on/off once every second.
13. If the battery voltage is below 6.6(+/-0.2) volts the unit will not turn on.
14. Open Electrode Detect: If an open circuit is detected at the output of channel A or B the output current will be reset at zero.

Expected service life:

5 years. Careful use and maintenance extends the life of the unit over the service life limit.

Physical dimensions: 108 x 65 x 23 mm

Weight: 70g without battery, 100g with battery.

Environmental Conditions for use:

+5 to +40 degrees Centigrade. 15-93% Humidity.

Environmental conditions for storage & transport:

-25 to +70 degrees Centigrade. 0-90% Humidity.



Issues and Solutions

Issue	Possible Cause	Solution
Battery symbol on LED flashing.	Low voltage.	Replace the battery.
Display does not come on and there is no signal from the unit.	Battery discharged.	Replace the battery.
	Battery incorrectly positioned.	Remove battery and reposition the battery in the compartment correctly.
	Battery contacts bent.	Use a tool (e.g., a standard small screwdriver) to push contacts in the battery compartment outward.
<p>Note: Depleted batteries are common, especially from above average use. Never tamper with the battery. Discard the battery if there is any indication of damage to the battery.</p> <p>Low Battery Indicator: If the battery's electromotive force goes below 6.9 (+/- 0.2) Volts the battery symbol will flash on/off once every second. If the battery voltage is below 6.6 (+/- 0.2) Volts the unit will not turn "ON". Dispose of batteries responsibly and in full compliance with all laws.</p>		
The unit turns "ON," but does not carry out commands (e.g., mA intensity increases but then drops to 0mA and readjustment to higher than 10mA level is not possible.).	Open circuit detected: Poor connection between the electrodes and your body.	Turn unit "OFF." Remove the PC Stim Electrode, and lubricate the electrode generously, reinsert the Electrode and then turn unit "ON."
	Open circuit detected: Lead wire connections are not secure.	Turn unit "OFF," unplug lead wires at all connection points and then replug all points, then turn unit "ON."
	Open circuit detected: Broken lead wires.	Turn unit "OFF," Replace lead wires or PC Stim Electrode unit, then turn unit "ON."

Issue	Possible Cause	Solution
Increasing intensity causes unpleasant sensation.	Dry spots or other factors that could increase resistance.	Turn unit “OFF.” Remove the PC Stim Electrode, lubricate it generously, reinsert and then turn unit “ON.”
		Use a different brand or style of lubricant
	Local muscle fatigue.	Be less aggressive in mA intensity or frequency of use.
<p>Never heat, freeze or microwave the PC Stim Electrode. Use only at eutheria (<i>e.g.</i>, normal healthy body temperature), 37.0° ± 0.5°degrees Celsius</p>		

Problem:

- Cannot reach maximum mA level; or
- The unit cuts off stimulation at certain level; or
- When increase the intensity, zero mA is flashing; or
- Power is cutting off when using

Solution:

It is normal behaviour in our and any other quality muscle stimulators (and TENS machines), and in most cases resolves itself - please read the guidance below.

The stimulation intensity will drop to zero if you simply press the mA+ button and no electrodes are connected to the channel on which you increase the intensity. You should attach a pair of electrodes to the lead wire and the lead should be connected to the channel on which you increase the stimulation intensity (mA).

Our unit is designed to detect any poor or intermittent connection across the electrodes and to cut off the stimulation output (mA) when it does so. This is a safety precaution. It is designed to prevent the user from inadvertently turning up the output stimulation current in the presence of a poor or intermittent connection and then experiencing a large unexpected powerful surge in the stimulation, if and when the connection is re-established.

Reasons for no connection if you use surface skin electrodes:

- * Check if both electrodes are connected to the same dual conductor lead wire, one electrode to the black connector (-) and another to red connector (+).
- * Be sure the electrodes adhere to your skin well and completely (e.g., not worn, unattached or curled at the edges). The electrode's conductive material is water-based and it will lose its adhesive qualities over time from perspiration, oil and debris the adhesive will collect. With care, maximum electrode lifespan can be achieved (see Care and Maintenance) but eventually the electrodes will need to be replaced.

At some point the electrodes will become dry. Moisten the adhesive surface with a few drops of water, and apply onto the plastic film overnight. This procedure will give you a few more days of electrode life.

If you find the electrodes will not stick due to oily skin, cleanse the skin with soap and water, then rinse and dry the area around the electrode site. If this does not work, try cleansing the skin with a swab impregnated with alcohol.

- * The most common reason for the device to detect a break in the circuit and return the mA setting to zero is a non-visible break in the lead wires (see the care and maintenance section of this manual). To check if the cable is good, cross the red and black pin and increase mA on the unit. If the cable conducts the electricity, the mA will go above 10mA and you would feel the stimulation mild tickling in your fingers which holds the crossed pins. If you feel a mild electrical current, this means the problem is with surface skin electrodes.



Warranty

DJO, LLC provides a warranty to the original purchaser, that this product will be free from defects in the material, components and workmanship, for a period of 2 years from the date of purchase by the distributor [invoice date from DJO to the appointed distributor]. If the distributor - from whom the product was purchased by the user - is satisfied that the product is defective, the user may return the unit directly to this distributor who will forward it to DJO. All such returns from the distributor to DJO must be authorised by DJO, in advance. The liability of DJO, under this limited product warranty does not extend to any misuse or abuse such as dropping or immersing the unit in water or other liquid substance or tampering with the unit or normal wear and tear. Any evidence of tampering will nullify this warranty.

Customer Service:

Please contact the selling dealer or DJO, LLC Service department.

Your invoice of purchase and/or the rear cover of this manual should state the name and the contact details of your distributor.

For assistance, if needed, in setting up, using or maintaining the unit, or report unexpected operation or events, please visit the manufacturer's website for further details: http://www.djoglobal.eu/fr_FR/index.html



MDSS GmbH, Schiffgraben 41
30175 Hannover, Germany

T: +49 511 6262 8630

F: +49 511 6262 8633



DJO, LLC

1430 Decision Street Vista, CA 92081 U.S.A

T: +1 760 727 1280

F: +1 760 734 5608

http://www.djoglobal.eu/fr_FR/index.html


This product is manufactured for DJO, LLC in compliance with the European Union Medical Device Directive MDD93/42/EEC under the supervision of BSI Group – EMEA, Notified Body number 0086.

DJO, LLC is certified by
BSI Group – EMEA to the
following Quality Standards: ISO13485:2003.




Dermatome Charts

Anterior View

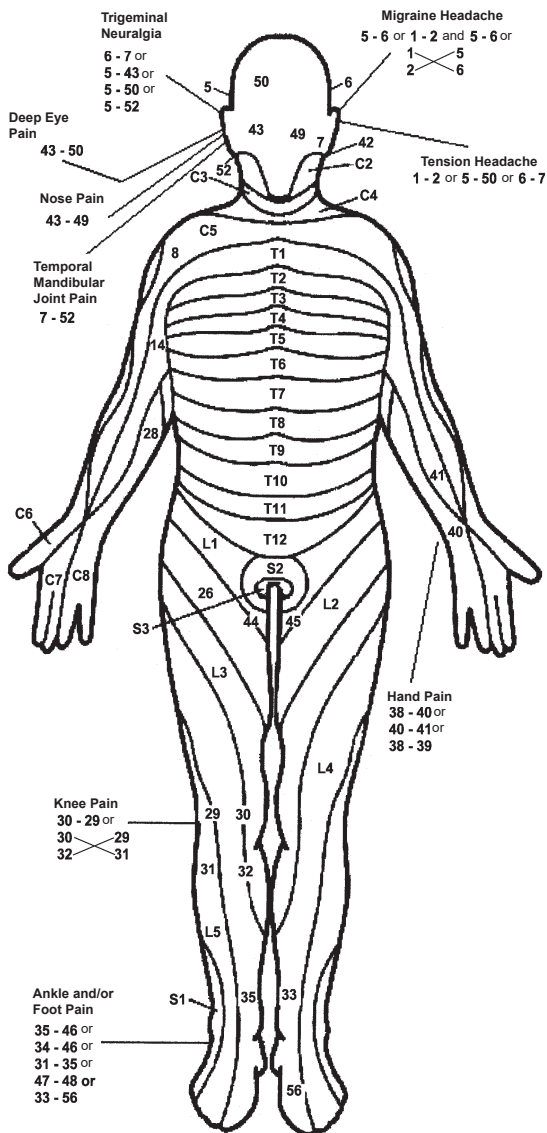

Posterior View.

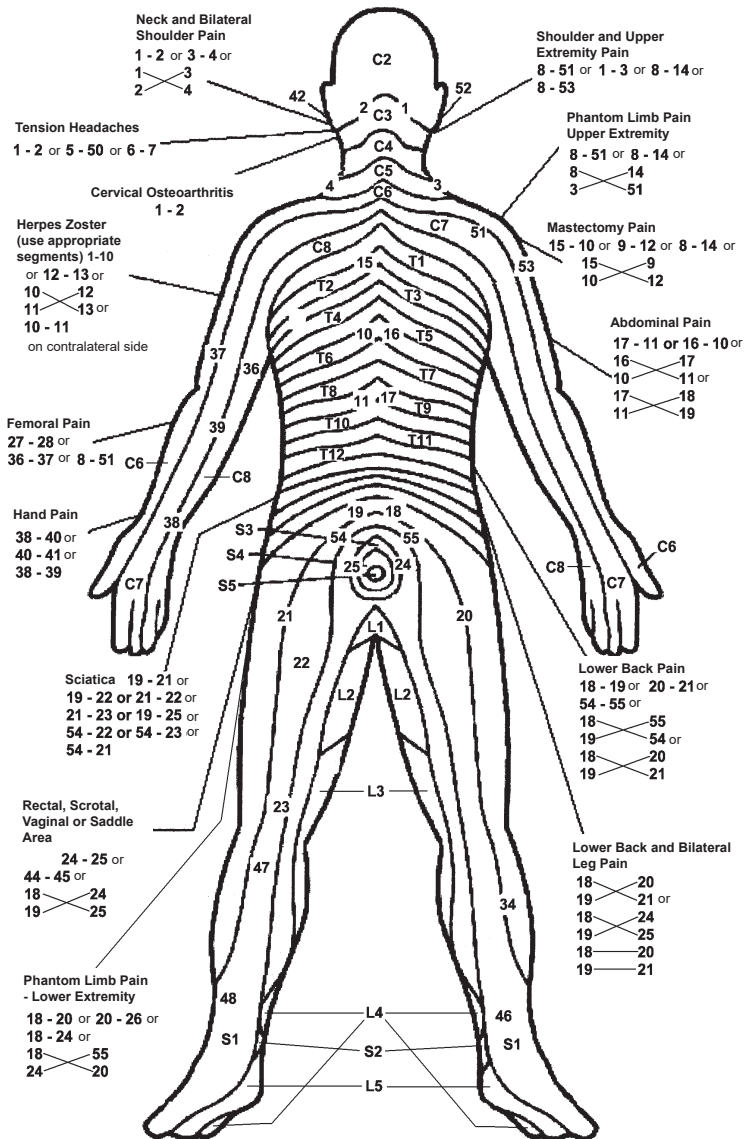

Clinical References

Please contact us for any clinical referencies of Primera, using the contact form on the website:

http://www.djoglobal.eu/fr_FR/index.html

User Profile

Patients, a patient's caretaker, or a family member providing assistance can use this device The user should be able to:

- Read and understand the directions, warnings and cautions
- Place the device on the patient
- Be able to see or hear device signals

Not for sale in the USA

Special Notes!

The Chattanooga Primera™ is equipped with a mA LOCKOUT.

Forty-five seconds after the unit is turned "ON", the intensity setting will lock. To make any adjustment to the intensity setting, you must first press the negative button to increase the mA.

Chattanooga Primera™



DJO, LLC
1430 Decision Street Vista,
CA 92081 U.S.A
T: +1 760 727 1280
F: +1 760 734 5608
http://www.djoglobal.eu/fr_FR/index.html

Document revision info.:

LOT 77622-OM-EN-K-08-02-16

Chattanooga
Primera International
Manual (English)

