





INTELECT® **RPW** LITE

EASY ACCESS TO RADIAL
PRESSURE WAVE THERAPY

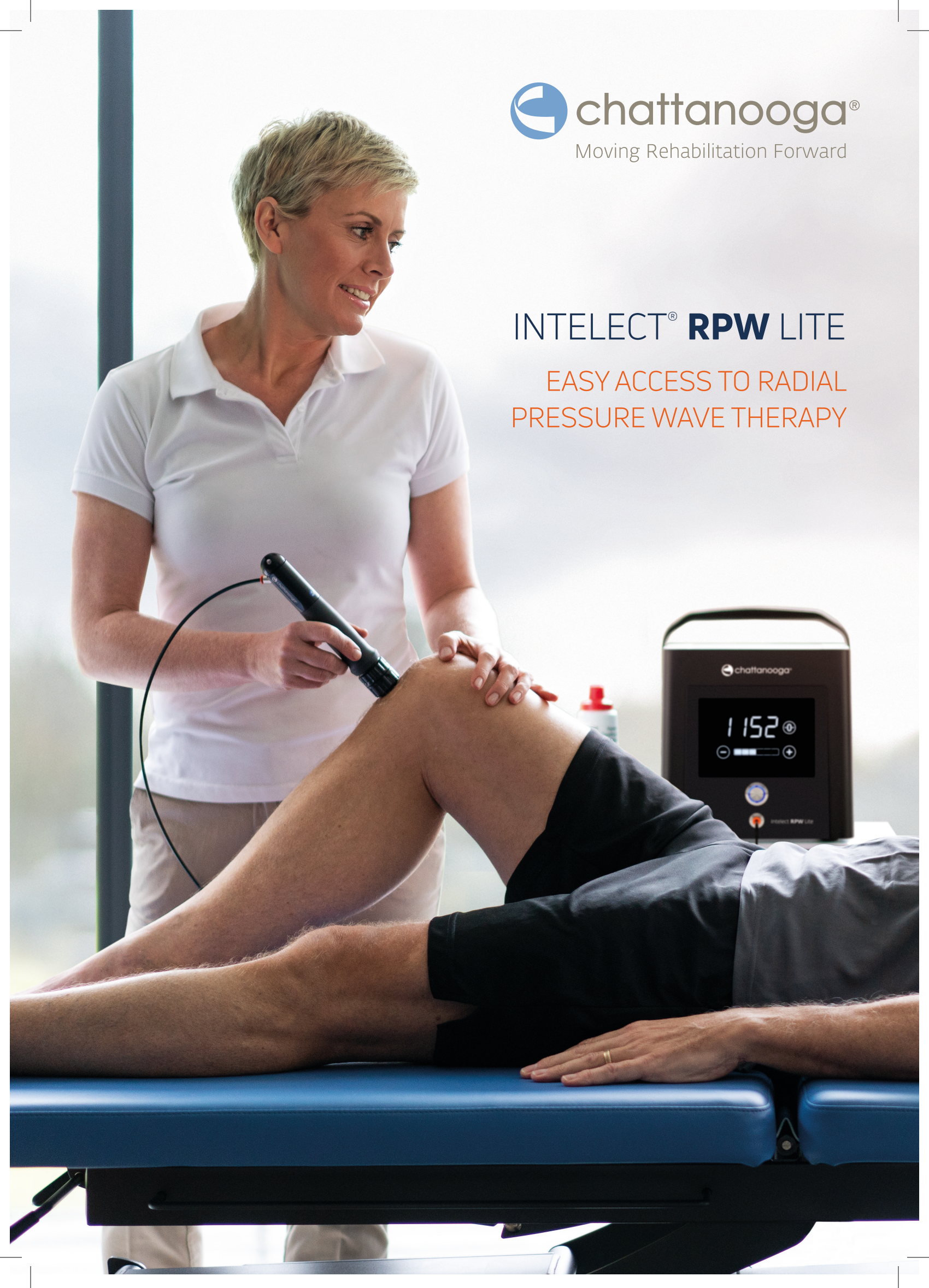

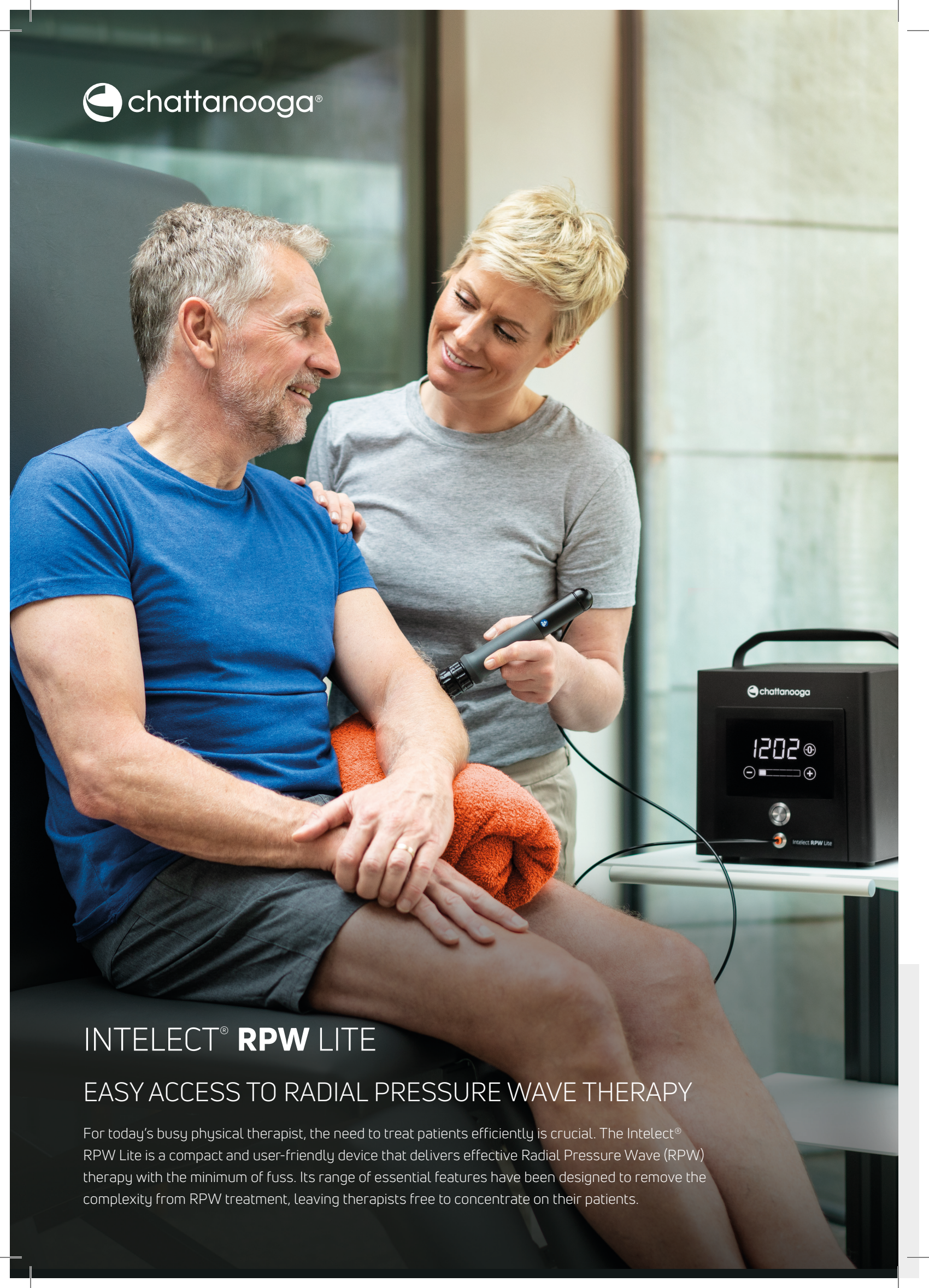

INTELECT® RPW LITE

EASY ACCESS TO RADIAL PRESSURE WAVE THERAPY

For today's busy physical therapist, the need to treat patients efficiently is crucial. The Intelect® RPW Lite is a compact and user-friendly device that delivers effective Radial Pressure Wave (RPW) therapy with the minimum of fuss. Its range of essential features have been designed to remove the complexity from RPW treatment, leaving therapists free to concentrate on their patients.

KEY FEATURES

PORTABILITY

Compact and easy to handle.

LARGE LCD TOUCH SCREEN

For quick and easy control of treatment parameters.

PULSE COUNTER

INTENSITY LEVEL


6 preset energy levels reduce set-up time whilst providing a useful range of settings for a variety of possible RPW indications.

EASY CONNECT

SPARROW™ handpiece technology with quick-release connection allows easy interchange of handpieces, ideal for treatments that require multiple different transmitters.



CHOICE OF 5 TRANSMITTERS

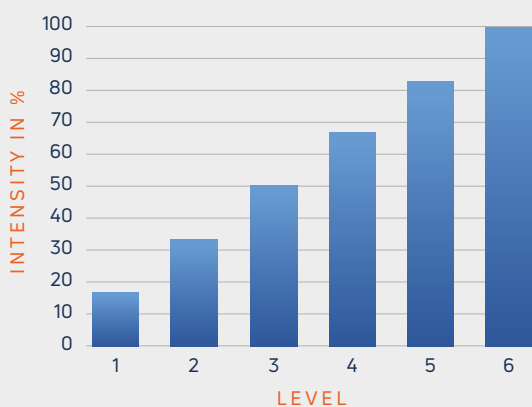

Relationship between intensity and frequency	
INTENSITY	FREQUENCY
1	18 Hz
2	16 Hz
3	10 Hz
4	9 Hz
5	8 Hz
6	6 Hz

The Intelect® RPW Lite offers an easy way to adjust the intensity level. Starting with a low level (step 1 at 18Hz) up to maximum energy release (step 6 at 6 Hz). The frequency changes automatically based on the intensity.

SPARROW™ is a registered trade mark of Storz medical AG

TARGETED PAIN RELIEF

The Intellect® RPW Lite allows for the treatment of indications with radial, pneumatically-generated, low-energy acoustic waves, or 'pressure pulses'. The applied pressure pulse propagates radially within the tissue, generating a therapeutic effect.

INDICATIONS

- ACHILLES TENDINOPATHY (1)
- CARPAL TUNNEL SYNDROME (2)
- PLANTAR FASCIITIS (3)
- TRIGGER FINGER (4)
- TRIGGER POINTS (5)
- PAIN AND INFLAMMATION MANAGEMENT IN ORTHOPAEDIC CONDITIONS (6)





CLINICAL INFORMATION

1. Gerdemeyer L. Current evidence of extracorporeal shock wave therapy in chronic Achilles tendinopathy. NCBINCB I Logo Int J Surg. 2015 Dec;24(Pt B):154-9. doi: 10.1016/j.ijssu.2015.07.718. Epub 2015 Aug 29.*

2. Wu YT. Effect of radial shock wave therapy for carpal tunnel syndrome: A prospective randomized, double-blind, placebo-controlled trial. J Orthop Res. 2016 Jun;34(6):977-84. doi: 10.1002/jor.23113. Epub 2015 Dec 10.

3. Ogden JA. Shockwave therapy for chronic proximal plantar fasciitis: a meta-analysis. Foot Ankle Int. 2002 Apr;23(4):301-8

4. Yildirim P. Extracorporeal shock wave therapy versus corticosteroid injection in the treatment of trigger finger: a randomized controlled study. J Hand Surg Eur Vol. 2016 Nov;41(9):977-983. doi: 10.1177/1753193415622733. Epub 2016 Sep 28.

5. Nedelka T. Mechano-transduction effect of shockwaves in the treatment of lumbar facet joint pain: comparative effectiveness evaluation of shockwave therapy, steroid injections and radiofrequency medial branch neurotomy. Neuro Endocrinol Lett. 2014;35(5):393-7.

6. Schmitz C. Efficacy and safety of extracorporeal shock wave therapy for orthopedic conditions: a systematic review on studies listed in the PEDro database. Br Med Bull. 2015;116:115-38. doi: 10.1093/bmb/ldv047. Epub 2015 Nov 18.

INTELECT® RPW LITE TRANSMITTERS

SPARROW™ transmitters are used with SPARROW handpiece only.

MYOFASCIAL TREATMENT



D20-S BLACK RPW LITE

Standard Oscillator, 20mm transmitter for muscle and connective tissue.

Penetration depth 0-50mm

Intensity level: Medium



D20-T BLACK RPW LITE

Golden Oscillator 20mm transmitter for professional athletes. High energy application.

Penetration depth 0-50mm

Intensity level: Very High

TRIGGER POINTS / TENDINOPATHIES



Ro40 BLACK RPW LITE

15 mm beam transmitter for any type of tendinopathy.

Penetration depth 0-35mm

Intensity level: Medium



C15 CERAmA-ax® BLACK RPW LITE

CERAmA-ax® Ceramic Energy 15mm transmitter for any type of tendonopathies.

Penetration depth 0-35mm

Intensity level: High



D115 BLACK RPW LITE

Golden Depth, 15mm Deep Impact® transmitter for deep target areas, chronic disorders and local trigger points.

Penetration depth 0-50mm

Intensity level: High



RADIAL PRESSURE WAVE TREATMENT

As a result of clinical evidence, most treatment protocols are based on 5 to 10 sessions of 1,500-2,500 shocks.

TREATMENT GUIDELINES

- Localize the painful areas you plan to treat
- Select the most appropriate transmitters based on treatment goal and targeted tissue
- Apply water-based ultrasound gel to the skin overlying the treatment area
- Always start with a lower level of intensity
- Perform a post-treatment test for improved function and decrease in pain



TECHNICAL INFORMATION

30680-230V: 230VAC 50-60Hz 700VA

30680-115V: 115VAC 60Hz 700VA

30680-100V: 100VAC 60Hz 700VA

DIMENSIONS: 289x238x310mm

WEIGHT: 9.8Kg

ENERGY: 6 Energy levels: 6-18Hz/max 2.7 bar

CONFORMITY: this device complies with the applicable standards EN/IEC 60601-1, CAN/CSA-C22.2 No.601.1, UL Std. No 60601-1.

CERTIFICATION: CE 0197

CLASSIFICATION ACCORDING TO MDD: Class IIa device 93/42 60601-1



INTELECT® RPW LITE

ORDERING INFORMATION

PART NUMBER	DESCRIPTION
30680-230V	Intelect® RPW Lite 230V set including standard accessories
30680-115V	Intelect RPW Lite 115V unit set including standard accessories
30680-100V	Intelect RPW Lite 100V unit set including standard accessories

STANDARD ACCESSORIES

PART NUMBER	DESCRIPTION
30885-1	Intelect RPW Lite SPARROW™ handpiece kit grey, including 1x projectile; 1x guiding tube; 1x wrench; 1x allen key
30873	Transmitter D20-S Black RPW Lite, 20mm
-	1x Power cord
30891	Operating manual and quick start guide
22601	Transmission gel, 250 ml
28739	Guiding tube cleaning brush
30890	Tube set (includes 3 tubes)

OPTIONAL ACCESSORIES

PART NUMBER	DESCRIPTION
30885-1	Intelect RPW Lite SPARROW handpiece kit grey, including 1x projectile; 1x guiding tube; 1x wrench; 1x allen key
30885-2	Intelect RPW Lite SPARROW handpiece kit green, including 1x projectile; 1x guiding tube; 1x wrench; 1x allen key
29104	SPARROW handpiece revision kit, including 2x projectiles; 2x guiding tubes short; 2x o-ring set; 1x o-ring guide
30873	Transmitter D20-S Black RPW Lite, 20mm
30876	Transmitter D20-T Black RPW Lite, 20mm
30842	Transmitter Ro40 Black RPW Lite, 15mm
30875	Transmitter C15 Black RPW Lite, 15mm
30874	Transmitter DI15 Black RPW Lite, Deep Impact 15mm
30883	Handpiece holder black
30890	Tube set (includes 3 tubes)
22601	Transmission gel, 250 ml
28739	Guiding tube cleaning brush
0.0032.012	Power cord EU
0.0032.012-US	Power cord USA
13-7458-CN	Power cord China/Australia/New Zealand
13448-CH	Power cord Switzerland



DJO® Global | 1a Guildford Business Park | Guildford | Surrey | GU2 8XG | UNITED KINGDOM

www.DJOGlobal.eu

SPARROW™ is a registered trade mark of Storz medical AG

Individual results may vary. Neither DJO Global, Inc. nor any of its subsidiaries dispense medical advice. The contents of this document do not constitute medical, legal, or any other type of professional advice. Information related to various health, medical, and fitness conditions and their treatment is not meant to be a substitute for the advice provided by a physician or other medical professional.

Copyright © [2019] by DJO, LLC All rights reserved. No part of this publication may be reproduced, distributed, or transmitted in any form or by any means.