



Hot and Cold

Introduction

Chattanooga is the world's largest manufacturer of rehabilitation equipment for treating musculoskeletal, neurological and soft tissue disorders. For over six decades, Chattanooga has set a benchmark for leadership, reliability and excellence. Our products contribute to better treatment outcomes in hospitals, clinics and home settings worldwide. We lead the physiotherapy industry by example and, through continuous innovation, provide real world solutions for clinicians and their patients.

As the world's leading brand in rehabilitation, Chattanooga is sold in over 80 countries worldwide. Chattanooga is directly represented in our key markets but also has long lasting relationships with dedicated distributors in the USA, Canada, Latin America, EMEA and Asia Pacific. Our Chattanooga distributors and customers are supported for order placement, training, product delivery and after sales service by a dedicated team of highly skilled DJO employees based all over the world.

Chattanooga is a brand of DJO Global Inc. DJO Global provides solutions for musculoskeletal health, vascular health and pain management. Our products help prevent injuries or rehabilitate after surgery, injury or from degenerative disease, enabling patients to regain or maintain their natural motion.

Visit www.DJOglobal.eu



Contents

Introduction	1
Hydrocollator Mobile Heating Units	3
Hydrocollator Stationary Heating Units	4
Hydrocollator HotPacs	5
Terry Covers	6
Stainless Steel Carts	6
ColPac Chilling Unit	7
ColPac, boo-boo pac, Sport-Pac, Nylatex Wraps	8



Hot and Cold

Chattanooga is renowned for its "hot and cold products" all over the world.

Chattanooga is the proud inventor of the original Hydrocollator HotPacs and heating unit, an all natural product that was first to market more than 60 years ago and is still the first choice of physiotherapists and other medical professionals when using heat therapy.

This HotPac generates deep moist heat which delivers increased advantages over dry heat to relieve chronic arthritic pain, acute muscle strains and general musculoskeletal discomfort.

As with moist heat therapy, Chattanooga sets the standard in cold therapy with the original ColPac range, which is now sold all over the world.

Ordering Information

Part Number Description

M-2 Mobile Hydrocollator

2402 120 V Unit

2407 220 V Unit

Standard Accessories

1006 Standard size HotPacs (12x)

M-4 Mobile Hydrocollator

2502 120 V Unit

2506 220 V Unit

Standard Accessories

1006 Standard size HotPacs (24x)

SS-2 Mobile Hydrocollator

2302 120 V Unit

2307 220 V Unit

Standard Accessories

1006 Standard size HotPacs (8x)

Technical Specifications

M-2 Mobile Hydrocollator

Mains Power: 110~240 V 50/60 Hz

Power Consumption: 1000 W

Weight: 39 kg

Dimensions*: 67 x 40 x 84 cm

Electrical Safety Class: Class 1, Type B

Safety Tests: EN 60601-1

M-4 Mobile Hydrocollator

Mains Power: 110~240 V 50/60 Hz

Power Consumption: 1500 W

Weight: 60 kg

Dimensions*: 89 x 51 x 84 cm

Electrical Safety Class: Class 1, Type B

Safety Tests: EN 60601-1

SS-2 Mobile Hydrocollator

Mains Power: 110~240 V 50/60 Hz

Power Consumption: 1000 W

Weight: 30 kg

Dimensions*: 53 x 41 x 84 cm

Electrical Safety Class: Class 1, Type B

Safety Tests: EN 60601-1

* Dimensions are listed as L x W x H

Hydrocollator® Mobile Heating Units



Features

Hydrocollator Mobile Heating Units

These durable stainless steel water tanks house the Hydrocollator Moist Heat HotPacs. The units are thermostatically controlled to ensure the ideal therapeutic temperature for HotPacs. Hydrocollator remains the benchmark heating unit against which all other units are judged.

- Full fiberglass insulation to prevent heat loss
- Dependable, rugged stainless steel design
- Simple to fill with water. No plumbing required
- Constant temperature of HotPacs is maintained
- Mobile units are equipped with 8 cm swivel, rubber casters for friction free movement about the clinic
- One year warranty

	M-2	M-4	SS-2	E-2	E-1
Tank Capacity	69 L	136 L	49 L	43 L	15 L
Temperature Range	71°-74°C	71°-74°C	71°-74°C	71°-74°C	71°-74°C
Thermal Cut-Out Temp.	82°-85°C	82°-85°C	82°-85°C	82°-85°C	82°-85°C
Accuracy	10%	10%	10%	10%	10%
Heat up Time (to 70°C)	8 Hours	8 Hours	4 Hours	4 Hours	2 Hours
Cool Down Time (from 70°C)	3 Hours	4 Hours	4 Hours	2 Hours	1 Hour
Fiberglass Insulation	Yes	Yes	Yes	No	No

Hydrocollator® Stationary Heating Units



E-1

E-2

Hydrocollator Accessories

Side Table Rack

Stainless steel work surface and wire shelves add efficiency to units. Racks available to fit most models (M-4, M-2).

Lotion Warmers

Keeps massage lotion and ultrasound gel at a comfortable temperature. Easily attaches to any treatment table or wall. Mini-Lotion Warmer hangs on the front of the Hydrocollator E-1 or E-2 heating units (bottles not included).

Wall Mounted Towel Rack

This stainless steel rack provides a convenient way to dry or store towels and terry covers (towels not included).



Ordering Information

Part Number Description

E-1 Stationary Hydrocollator

2102	120 V Unit
2107	220 V Unit

Standard Accessories

1006	Standard size HotPacs (4x)
------	----------------------------

E-2 Stationary Hydrocollator

2802	120 V Unit
2806	220 V Unit

Standard Accessories

1006	Standard size HotPacs (3x)
1004	Oversize HotPacs (2x)
1002	Neck Contour HotPac (1x)

Hydrocollator Accessories

Side Table Rack

4010	Side Table Rack (not shown)
4230	Extra Shelf for Side Table Rack

TM-1 Electric Lotion Warmer

4012	TM-1 Electric Lotion Warmer (110 V only)
------	--

Mini Non-Electric Lotion Warmer

22100	Mini Lotion Warmer
-------	--------------------

Dial Thermometer

4228	Dial Thermometer (not shown)
------	------------------------------

Wall Mounted Towel Rack

4016	Wall Mounted Towel Rack 5 x 84 cm
------	-----------------------------------

Technical Specifications

E-1 Stationary Hydrocollator

Mains Power:	110~240 V 50/60 Hz
Power Consumption:	1000 W
Weight:	7 kg
Dimensions*:	33 x 20 x 41 cm
Electrical Safety Class:	Class 1, Type B
Safety Tests:	EN 60601-1

E-2 Stationary Hydrocollator

Mains Power:	110~240 V 50/60 Hz
Power Consumption:	1000 W
Weight:	14 kg
Dimensions*:	38 x 33 x 51 cm
Electrical Safety Class:	Class 1, Type B
Safety Tests:	EN 60601-1

* Dimensions are listed as L x W x H


Ordering Information

Part Number Description

Hydrocollator HotPac - Moist Heat Hot Pack

A 1008	23 cm x 39 cm
B 1010	25 cm x 61 cm
C 1004	Oversize 38 cm x 61 cm
D 1002	Neck Contour 61 cm
E 1014	Half Size 13 cm x 30 cm
F 1006	Standard 25 cm x 30 cm
G 1012	Knee-Shoulder 25 cm x 50 cm

Hydrocollator® HotPacs™

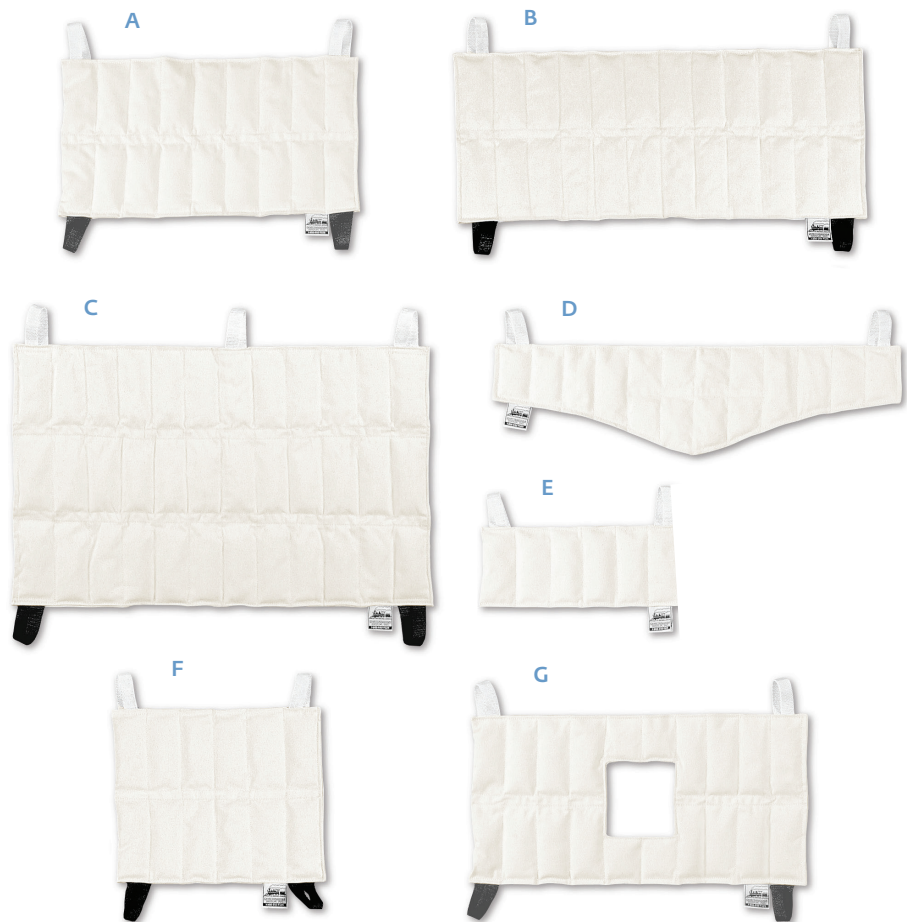

Features

The Original Hydrocollator HotPac

Hydrocollator Moist Heat HotPacs are the first choice for physiotherapists, physicians, chiropractors, athletic trainers and kinesiologists for over 60 years. Chattanooga is the inventor of the original Hydrocollator HotPac and proud to still be producing this all natural product today. The benefits of moist heat therapy over dry heat are well established for pain relief due to arthritis, myofibrosis, back pain, muscle strains and spasms.

Hydrocollator HotPacs

- Seven HotPac sizes designed to fit every body contour
- Provides up to 30 minutes of deep, soothing, moist heat
- Tabs on most HotPacs for easy rotation in busy clinics
- Cost effective, provides hundreds of treatments from one HotPac
- All natural pack containing bentonite, a natural hydrocolloid filling, wrapped in a woven cotton blend
- One year warranty



Terry Covers



Features

All-Terry Covers

- For high-volume institutional use
- Hook and loop closures to prevent unwrapping or slipping

Foam-Filled Terry Covers

- For home or clinic use
- Polyurethane foam layer between two layers of terry
- Available with or without pockets

Ordering Information

Part Number Description

All-Terry Covers

1118	Standard 48 cm x 69 cm
1126	Oversize 61 cm x 91 cm
1104	Neck Contour 64 cm x 46 cm

Foam-Filled Terry Covers

1102	Standard Size 48 cm x 69 cm
1124	Oversize 61 cm x 91 cm
1120	Neck Contour 64 cm x 41 cm

Foam-Filled Pocket Terry Covers

1108	Standard 43 cm x 69 cm
1110	Oversize Pocket 61 cm x 91 cm

Hydrocollator HotPac Sets (Combines HotPac and Foam Filled Terry Cover with retail packaging)

1064	Standard Set
1062	Neck Set

Stainless Steel Carts

4002	Model MB, Specially designed to hold the E-2 and E-1 Hydrocollator heating units 51 x 48 x 43 cm (20"H x 19"W x 17"D)
4004	Model MB-T, Accommodates stationary heating, ultrasound and electrotherapy units 76 x 48 x 43 cm (30"H x 19"W x 17"D)
4018	Model MB-TD, Accommodates stationary heating and electrotherapy units. Supply drawer included 76 x 48 x 43 cm (30"H x 19"W x 17"D)

Stainless Steel Carts



Ordering Information

Part Number Description

ColPac Chilling Unit Model C-5

3102 110 V Unit

3114 220 V Unit

Standard Accessories

1500 Standard Size ColPacs (6x)

1506 Half Size ColPacs (6x)

Technical Specifications

ColPac Chilling Unit

Mains Power: 220~240 V 50 Hz; 115 V, 60 Hz

Weight: 74 kg

Dimensions*: 69 x 41 x 84 cm

Electrical Safety Class: Type B

Safety Tests: EN 60601-1

* Dimensions are listed as L x W x H

ColPac® Chilling Unit



Features

ColPac Chilling Unit

Large capacity and faster chilling time mean your patients will never have to wait for a ColPac. The ColPac Chilling Unit helps you keep pace throughout a busy day. No plumbing required, simply plug into a standard electrical outlet.

All ColPac Chilling Units include:

- Closed-cell foam insulation for efficiency
- Stainless steel, welded construction for long life
- Heavy-duty compressor for more efficient cooling
- Specially designed interior chilling coils cool packs faster
- 8 cm (3") swivel-type rubber casters for silent, friction-free movement of mobile units
- Drain valve for easy cleaning and defrosting
- One year warranty



Colpac®, boo-boo pac™, Sport-Pac™, Nylatex® Wraps



Features

The Original ColPac

Cold therapy, when applied directly to the affected area, is used to control pain due to inflammation, swelling, multiple sclerosis, muscle strains, sprains and spasms.

- Seven different sizes of blue vinyl ensure effective cold therapy regardless of body area or shape. Also available are four sizes of black polyurethane
- Effective relief of pain
- Filled with a non-toxic silica gel that remains pliable throughout treatment
- Delivers up to 30 minutes of soothing relief from acute pain, swelling and fever
- Latex free
- One year warranty

boo-boo pac

- Fun, bear-shaped design especially made for kids
- Ideal for pediatric applications
- Filled with non-toxic silica gel
- Delivers up to 30 minutes of soothing relief from acute pain, swelling and fever

Nylatex Wraps

- Ideal for holding ColPacs in place
- Several sizes available
- Good for hundreds of uses

FlexiPac

- Microwaveable to provide soothing heat, or after freezer storage to provide penetrating cold therapy
- Flexible gel conforms to any body contour at any temperature

Sport-Pac

- Soccer ball design is ideal for youth sport applications
- Easily transferred from the freezer to the team cooler for sideline use
- Provides 30 minutes of effective cold therapy for acute pain due to strains, sprains and bruises

Ordering Information

Part Number Description

Blue Vinyl ColPac

A 1512	Oversize 28 cm x 53 cm
B 1500	Standard Size 28 cm x 36 cm
C 1506	Half Size 19 cm x 28 cm
D 1508	Neck contour 58 cm
E 1502	8 cm x 28 cm
F 1504	Quarter Size 14 cm x 19 cm
G 1510	Eye (Specially shaped to cover closed eyelids)

Black Polyurethane ColPac

H 1554	Neck Contour 53 cm Long
I 1556	Oversize 32 cm x 47 cm
J 1552	Standard Size 25 cm x 34 cm
K 1562	Half Size 17 cm x 28 cm

L boo-boo pac (Shipped in a display box containing 12 packs)

1534	Royal Blue boo-boo pac
------	------------------------

M Nylatex Assortment (1 pack of 9)

1200	Nylatex Assortment (One of each size listed below except 6 cm x 61 cm)
------	--

Nylatex Individual (1 pack of 3)

1202	6 x 46 cm
1204	6 x 91 cm
1206	6 x 122 cm
1208	10 x 46 cm
1210	10 x 91 cm
1212	10 x 122 cm
1213	10 x 152 cm
1214	15 x 46 cm
1216	15 x 91 cm
1218	15 x 122 cm
1219	15 x 152 cm
1264	6 x 61 cm

N FlexiPac Hot & Cold Compresses

4026	FlexiPac 3 cm x 15 cm (48/case)
4020	FlexiPac 13 cm x 25 cm (24/case)
4029	FlexiPac 20 cm x 35 cm (12/case)

O Sport-Pac (Shipped in a display box containing 10 packs)

1517	Sport-Pac
------	-----------

DJO GLOBAL

AUSTRALIA:

T: +1300 66 77 30

F: +1300 66 77 40

E: customerservice.au@DJOglobal.com

BENELUX:

T: Belgium 0800 18 246

T: Netherlands 0800 0229442

T: Luxemburg 8002 27 42

E: benelux.orders@DJOglobal.com

CANADA:

T: +1 866 866 5031

F: +1 866 866 5032

E: canada.orders@DJOglobal.com

CHINA:

T: (8621) 6031 9989

F: (8621) 6031 9709

E: information_china@DJOglobal.com

DENMARK, FINLAND,

NORWAY & SWEDEN:

T: Sweden 040 39 40 00

T: Norway 8006 1052

T: Finland 0800 114 582

T: Denmark +46 40 39 40 00

E: info.nordic@DJOglobal.com

FRANCE:

T: +33 (0)5 59 52 86 90

F: +33 (0)5 59 52 86 91

E: sce.cial@DJOglobal.com

GERMANY:

T: +49 761 4566 01

F: +49 761 456655 01

E: infoservice@DJOglobal.com

ITALY:

T: +39 02 484 63386

F: +39 02 484 09217

E: vendite@DJOglobal.com

INDIA:

T: +91 44 66915127

E: customercare.india@DJOglobal.com

SOUTH AFRICA:

T: +27 (0) 87 3102480

F: +27 (0) 86 6098891

E: chattanooga.saf@DJOglobal.com

SPAIN:

T: +34 934 803 202

F: +34 934 733 667

E: ventas@DJOglobal.com

SWITZERLAND:

T: +41 (0) 21 695 2360

F: +41 (0) 21 695 2361

E: info@compex.ch

UK & IRELAND:

T: +44 (0)1483 459 659

F: +44 (0)1483 459 470

E: ukorders@DJOglobal.com

UNITED STATES:

T: +1 800 336 6569

F: +1 800 936 6569

E: customercare@DJOglobal.com

DJO GLOBAL, EXPORT CENTRES

ASIA-PACIFIC:

DJO Asia-Pacific Limited

Unit 1905, 19/F, Tower II

Grand Central Plaza

138 Shatin Rural Committee Road

Shatin

HONG KONG

T: +852 3105 2237

F: +852 3105 1444

E: info.asia@DJOglobal.com

EUROPE, MIDDLE EAST & AFRICA:

DJO Benelux

Welvaartstraat 8

2200 Herentals

BELGIUM

T: +32 (0) 14248350

F: +32 (0) 14248358

E: info.emea@DJOglobal.com

LATIN AMERICA:

DJO Global, Inc

1430 Decision Street

Vista

CA 92081-8553

U.S.A.

T: 1 800 336 6569

F: 1 800 936 6569

E: info.latam@DJOglobal.com



⁺
MOTION IS MEDICINE™