

Ordering Information	
Part Number	Description
Montane Tattras - 46cm Head Section	
3520	Tattras 2 Section 120 V - 2 Wheels / 2 Legs
3521	Tattras 2 Section 120 V - 4 Wheels
3522	Tattras 2 Section 230 V - 2 Wheels / 2 Legs
3523	Tattras 2 Section 230 V - 4 Wheels
Montane Tattras Gerlach - 69cm Head Section	
3525	Montane Tattras Gerlach 2 Section 230 V - 4 Wheels

Technical Specifications	
Mains Power:	230 V / 50 Hz, 120 V / 60 Hz
Standby Consumption:	35 W
Input:	230 V / 1.0 A / 50 Hz, 120 V / 1.9 A / 60 Hz
Weight:	110 kg (242 lbs) plus 30 kg pallet when shipping
Dimensions*:	198 x 69 x 49 - 101 cm
Shipping Dimensions*:	219 x 84 x 75 cm
Lifting Capacity:	150 kg (330 lbs)
Electrical Safety Class:	Class II, Type B
Safety Tests:	Meets Directive 93 / 42 / EEC, Certified to CAN/CSA Standard C22.2 No. 601.1-M90x/A2
* Dimensions are listed as L x W x H	

Upholstery Color Options		
		
Ecru-Crème (x = 00)	Beige (x = 01)	Blue (x = 02)
		
Navy Blue (x = 03)	Terra Cotta (x = 04)	Chocolate (x = 05)
		
Gray (x = 06)	Graphite Gray (x = 07)	Black (x = 08)

NOTICE: Color charts are produced within the limitations of lithographic color printing and are intended as a general guide only. For ordering purposes, please use model number, hyphen, upholstery color number as indicated in the chart above. Example: Montane Andes 7 Postural Drainage in Chocolate - 3578-05

Montane Tattras/Tattras Gerlach 2 Section Table



Montane Tattras



Montane Tattras Gerlach




Features

- This table is available in 2 configurations - A 46cm head section and a 69cm head section option to satisfy different clinical needs
- Electric high/low height adjustment from 49 to 101 cm for easy and safe patient transfer and adjustment to the ideal working height for treatment
- Easy access foot switch, accessible from any position around the table 360°, for simple height adjustment
- Features smooth vertical lift with a capacity up to 150 kg
- Adjustable head support - positive/negative inclinations
- Face hole on head section with removable cushion
- Side brackets for stabilizing belts are standard equipment

Table Articulations

Montane Tattras

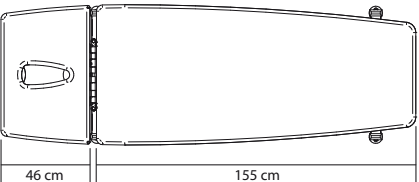

Montane Tattras Gerlach

